

 **smeg**
INSTRUMENTS



SMEG PROFESSIONAL INSTRUMENT WASHERS

Accessories

The GW3050H, 4050H and 5050H machines must be equipped with the most suitable trolleys and supports for the different types of surgical instrument used. The following descriptions of the different types of machine/trolley configurations can serve as a guide for choosing the most suitable washing solution.

SURGICAL INSTRUMENTS WASHING



CSK4

- 2 WASHING LEVELS

CSK4 SURGICAL INSTRUMENT TROLLEY

Washing trolley for surgical instruments. Suitable for treating standard and large-sized surgical instruments. Two washing levels served by an incorporated sprayer and a sprayer on the bottom of the machine. To be fitted with N.4 CSK2 baskets. Made of stainless steel.

CSK2 SURGICAL INSTRUMENT BASKET

Surgical instrument basket with handles. Occupies 1/2 level of CSK4 and CSK8. Dimension LxDxH 450x225x50 mm., mesh 5x5 mm. Made of stainless steel.



CSK6

- 3 WASHING LEVELS

CSK6 SURGICAL INSTRUMENT TROLLEY

Washing trolley for surgical instruments. Suitable for treating large quantities of standard-sized surgical instruments. Three washing levels served by two incorporated sprayers and one sprayer on the bottom of the machine. To be fitted with 6 CSKDIN baskets. Made of stainless steel.

CSKDIN SURGICAL INSTRUMENT BASKET

Surgical instruments basket with handles. Occupies 1/2 level of CSK6. Dimension LxDxH 480x245x50 mm., mesh 5x5 mm. Made of stainless steel.

- 4 WASHING LEVELS

CSK8 SURGICAL INSTRUMENT TROLLEY

Washing trolley for surgical instruments. Suitable for treating standard and large-sized surgical instruments. Four washing levels served by three incorporated sprayers and one sprayer on the bottom of the machine. To be fitted with 8 CSK2 baskets. The last solution is particularly suitable for treating large quantities of instruments from a variety of sources, which have to be kept separate. Made of stainless steel.



CSK8

CSK2 SURGICAL INSTRUMENT BASKET

Surgical instrument basket with handles. Occupies 1/2 level of CSK4 and CSK8. Dimension LxDxH 450x225x50 mm., mesh 5x5 mm. Made of stainless steel.

CSK1 SURGICAL INSTRUMENT BASKET

Surgical instruments basket with handles. Suitable for treating large-sized surgical instruments (orthopaedic). Occupies 1 complete level of CSK4, CSK6 and CSK8. Dimension LxDxH 450x450x55 mm., mesh 5x5 mm. Made of stainless steel.

CSMICRO1 SURGICAL INSTRUMENT BASKET

Surgical instruments basket with handles and cover. Suitable for treating tiny surgical instruments. Occupies 1/4 level of CSK4, CSK6 and CSK8. Dimension LxDxH 240x235x55 mm., mesh 2x2 mm. Made of stainless steel.

RIGID ENDOSCOPE WASHING

EF2 RIGID ENDOSCOPE TROLLEY

Injection trolley for washing and drying of instruments for minimal invasive surgery. Provided with 38 injectors of various diameters, 10 special connectors and 3 accessory basket (n. 2 CSMICRO1 + n. 1 CSK2). Suitable for treating hollow instruments to be fitted to the nozzles or special connectors which inject the detergent and disinfectant solution into the channels inside the instruments. Connector with Drying System coupling valve for drying the inside of the instruments. Made of stainless steel.



EF2

Accessories



OPHTHALMOLOGY WASHING

EF3 TROLLEY FOR OPHTHALMOLOGY

Capacity: 21 positions for washing of hollow instruments, 6 catheters, 2 positions for external washing, 4 baskets positioned on two levels for washing various types of instruments by means of the spray arm. Connector with Drying System coupling valve for drying the inside of the instruments.

Made of stainless steel.

MINIMAL INVASIVE SURGERY WASHING

EF4 TROLLEY FOR MINIMAL INVASIVE SURGERY

Trolley for washing the internal parts of hollow instruments with 40 connections placed on rows: two rows for instruments placed vertically, two rows for instruments placed obliquely. Completed with spigots, supports and adapters of various dimensions for washing of instruments that can be dismantled, rigid optics etc; 1 CSK2 basket + 3 CSK1/6 BASKET included. Connector with Drying System coupling valve for drying the inside of the instruments. Made of stainless steel



ZCS3

OPERATING THEATRE SHOE WASHING

ZCS3 OPERATING THEATRE SHOE TROLLEY

Trolley with two washing levels for treating operating theatre shoes. It can hold 12 pairs. If other instruments, besides the shoes, are treated in the same instrument washer (using the special baskets/supports), the ZCS3 must be used with the Kit Z5000. Made of stainless steel.

Z5000 SPARE SPRAYER AND FILTER KIT

Bottom sprayer and filters to be fitted with trolley ZCS3 in multipurpose instrument washers.




ZCS3

ANAESTHESIA WASHING

AE2 ANAESTHESIA TROLLEY

Internal injection trolley for treating anaesthesia instruments. Capacity 8 flexible tubes of max. length 150 mm., 13 catheters, 3 breathing bags, 4 masks and accessories. Connector with Drying System coupling valve for drying the inside of the tubes. Made of stainless steel.

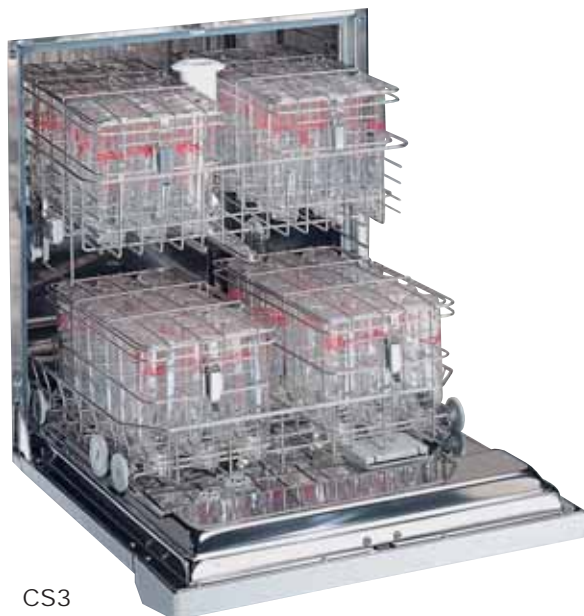

CS3

BABY BOTTLES WASHING

CS3 STANDARD TROLLEY

Standard trolley on two washing levels served by an incorporated sprayer and the machine bottom sprayer. Suitable for positioning of various supports.

Set for treating baby bottles: 4 CB (18-18/100-18/120)

Set for treating baby bottles and teats: 1 CPS20 + 3 CB (18-18/100-18/120). Made of stainless steel.

CB18/100 CONTAINER FOR 100 ML. BABY BOTTLES

Container for 18 x 100-ml. baby bottles with lid. Positioning on CS3 standard trolley. Occupies 1/2 washing level. Made of stainless steel.

CB18/120 CONTAINER FOR 120 ML. BABY BOTTLES

Container for 18 x 120-ml. baby bottles with lid. Positioning on CS3 standard trolley. Occupies 1/2 washing level. Made of stainless steel.

CB18 CONTAINER FOR 250 ML. BABY BOTTLES

Container for 18 x 250-ml. baby bottles with lid. Positioning on CS3 standard trolley. Occupies 1/2 washing level. Made of stainless steel.

CPS20 CONTAINER FOR 20 DUMMIES

Container for 20 dummies. Positioning on CS3 standard trolley. Occupies 1/2 washing level. Made of stainless steel.

SA3 KIDNEY BOWL SUPPORT

Basin and kidney bowl support. Positioning on CS3 standard trolley. Occupies 1/2 washing level. Made of stainless steel.

VARIOUS INSTRUMENTS WASHING

LBT5 BOTTLE TROLLEY

Injection trolley for washing 5 secretion bottles. Fitted with special nozzles for internal spraying.

Made of stainless steel.

LBT5DS

Version fitted with connector with Drying System coupling valve for drying inside bottles. Made of stainless steel.

CPC1 CARRELLO CONTAINERS

Trolley for washing 2 containers of dimensions LxDxH 300x300x150 mm. + lids. Two washing levels served by sprayers. Made of stainless steel.